ITFIT MOBILE LIFESTYLE

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備註: 1."RoHS 1.0" and "RoHS 2.0" indicate that the percentage content of the restricted substance(s) listed in the corresponding column of the table is less than 0.1% by weight of homogeneous material. 2."RoHS 2.0" indicates that the restricted substance(s) listed in the corresponding column of the table is less than 0.01% by weight of homogeneous material. 3."RoHS 2.0" indicates that the restricted substance(s) listed in the corresponding column of the table is less than 0.001% by weight of homogeneous material.

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ITFIT MOBILE LIFESTYLE

Before Use

Thank you for using Rechargeable Electric Toothbrush. Before operating this product, please read the User's Manual completely and save them for future use.

Feature of Product

- 5 Cleaning Modes
- 2-Minute Auto Shut-off
- DuPont Soft Brush Heads
- Rechargeable Battery
- IPX7 Waterproof Toothbrush
- Low-Noise Motor

Product Information


Product Name:ITFIT Rechargeable Electric Toothbrush
 Model No. ITFITEX38
 Material:ABS, Cylindrical Lithium-ion Battery, Electronic Parts
 Input(Charging Station):5V/1A
 Power Consumption(Charging Station):1.8W
 Battery Capacity:800mAh(2.96Wh)/3.7V
 Product Dimensions(Brush Head and Handle):245x30x29mm
 Product Weight(Brush Head and Handle):Approx. 116g
 Travel Case Dimensions:215x73x40mm
 Travel Case Weight:108g
 Operating/Storage Temperature:0°C-45°C
 Waterproof Rating:IPX7
 Country of Manufacture:China
 Manufacturing Date:2024.JUN

Rechargeable Electric Toothbrush User Manual

Contains

Toothbrush Handle, Brush Heads x 6pcs(Random Colors), Charging Station, Travel Case, User Manual

Product Appearance Introduction

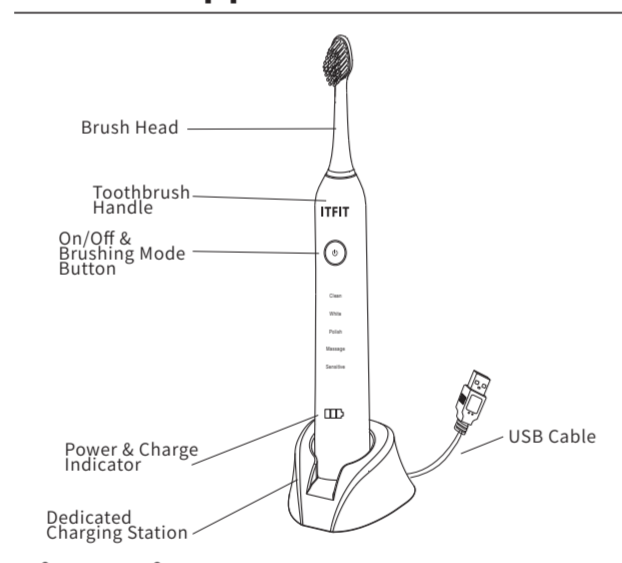• Mode Descriptions

Mode Names	Indicator Light	VPM(Vibrations Per Minute)	Vibration Frequency
Clean	Clean	31000	258±18Hz
White	White	41000	342±24Hz
Polish	Polish	48000	400±28Hz
Massage	Massage	31000	258±18Hz
Sensitive	Sensitive	31000	258±18Hz

• Briefly press the power button for manual shutdown (automatic shutdown after 2 minutes).
 Notes:
 • Within the initial 5 seconds before the mode starts, a short press of the power button shifts to the next mode. Beyond 5 seconds into the mode, a short press of the power button powers off the device.
 • Clean/Message/Sensitive modes share the same VPM and frequency but vary slightly in vibration intensity and pattern.
 • The toothbrush features a built-in timer with a prompt every 30 seconds, reminding you to clean each quadrant of your mouth for approximately a quarter of the recommended brushing time.

• First Time Use

- Store the brush head on the top of the handle after it is fully charged (maintain approximately 1mm gap between the brush head and the handle).
- Charge the product before initial use (requires 24 hours for a full charge from zero).
- Connect the charging station to a USB-A port for charging.
- Place the handle vertically on the charging station.
- Disconnect the power cord after the charge is complete.

• How to Use

- Dampen the brush head and apply toothpaste.
- To prevent toothpaste splatter, place the brush head inside the mouth before turning on the power.
- Press the power button briefly to start.
- After powering on, short press the power button to switch modes.

• Electric Toothbrush Usage Tips

- Start brushing from the edges of the teeth, gently guiding the electric toothbrush across the entire tooth.
- Avoid applying excessive pressure with the brush head to reduce the risk of tooth sensitivity.
- It's recommended to brush the lower teeth first, then the upper teeth. Brush the outer surfaces of the teeth before moving to the inner surfaces.
- Brush 2-3 teeth at a time during the brushing process.
- For a better use of the product, you may consult with a dentist.

• Charging Indicator Light Guide

Placing the handle on the charging station initiates the charging process, and the battery indicator lights reflect the battery status as follows:

Daily Maintenance

- Keep the brush head clean to prevent the growth of bacteria and microbes.
- After each use, rinse the brush head with water and pat it dry.
- It is recommended to detach the brush head from the handle and rinse it at least once a week.
- Regularly replace the brush head (recommended every 3 months). If you notice deformities or wear on the bristles, replace them promptly.
- Rinse the brush head and handle with water. Use a soft cloth to wipe the charging station and travel case.

Eco-friendly Tips

This product contains a battery and is considered recyclable electronic waste. To protect the environment, please do not dispose of it with household garbage; instead, take it to the electronic waste recycling facility provided in your area.

Warnings

- Consult a dentist before use if undergoing recent oral treatment.
- Ensure the brush head is securely attached to the handle before each use; do not use without the brush head.
- Do not use the product if unfamiliar with its operation to prevent potential harm.

Keep out of reach of children; not intended for use as a toy.

• Cease use immediately if the product, charger, or charging device is damaged.

• Avoid submerging the charger in water and refrain from prolonged exposure of the electric toothbrush to water to prevent malfunction, fire, electric shock, or explosion.

• If the charger accidentally falls into water, unplug it before retrieving.

• Do not unplug the device with wet hands.

• Avoid storing the product in high-temperature or humid areas.

• Do not clean the product in a dishwasher.

For brush head purchases, please contact our company.

Caution

- In the initial days of product use, slight gum bleeding may occur, which typically subsides within a week. If bleeding persists or any discomfort arises after a week, discontinue use and consult a dentist immediately.
- Do not dismantle or alter any parts of the product to prevent damage.
- This product is designed for personal oral care only; do not use it for other purposes.
- Not suitable for children, infants, or individuals unable to operate it unless under the supervision or guidance of a guardian.
- Do not place coins, clips, or other metal objects on the charging station.

Warranty Instructions

Warranty Period Description:
 For VIP gifts such as mobile phones/tablets sold by Samsung Electronics Taiwan, you can enjoy a 1-year warranty from the date of purchase with this warranty

and the purchase invoice of Samsung products.
 The free service will not be provided in the following cases:
 1.The product has exceeded the warranty period provided by the company.
 2.Damage or failure caused by natural disasters, insect (rat)damage, lightning strikes and other factors.
 3.Failures caused by other products (equipment) other than the machine itself.
 4.Failure to operate in accordance with product specifications and user manuals, resulting in product abnormalities or operating system changes.
 5.The cause of the failure is beyond the scope of the original warranty, such as: damp, damaged, deformed, using a non-original authorized repair center for the product or using non-original spare parts.
 6.Human failure during movement, transportation or improper placement.
 7.Failures caused by using other equipment and accessories that do not meet the applicable specifications of the product.
 8.Unsuitable power circuit and voltage and current.
 9.Appearance and consumable materials.
Declaration of Restricted Substances Content

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