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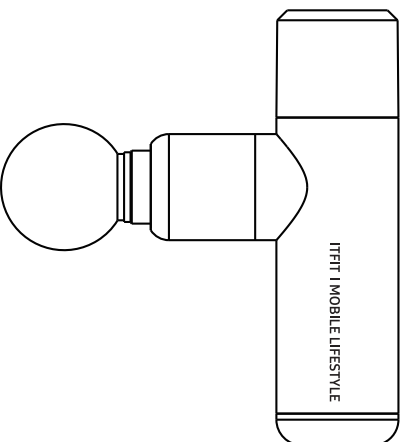
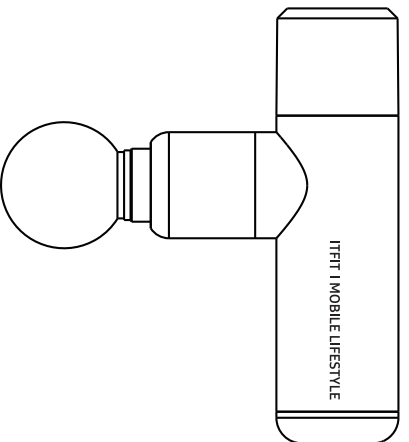
**Massage Gun
User Manual**

Before Use

Before operating this product, please read the User's Manual completely and save them for future use. Please keep the User's Manual in a safe place for reference at any time.

Feature of Product

- 4-speed massage with variable settings
- Can be used for up to 7 hours
- 10-minute automatic shutdown safety setting
- Durable metal body
- Compact and portable design
- Ergonomic handle for a comfortable grip


Product Information

Product Name: Massage Gun
Model No. ITFITEX43A
Materials: Metal, Electronic Parts
Battery Type: Li-ion Battery
Battery Capacity: 1200mAh 7.4V(8.88Wh)
Charging Time: Approx. 2.5 hrs.
Battery Life: Up To Approx. 7 hrs. of Continuous Use (Varies Depending on The Speed Setting)
Charging Type: USB Type C
Input(charge): 5V=2A
Speed Settings: 4 Speeds,
Ranging From 700 to 3000 Revolutions per Min.
Dimensions: 115 x 131 x 37mm
Net Weight: 311g
Production Date: 2024.11
Made in China

Contains

Pocket-sized Massage Gun, Massage Heads x 4, USB-C Charging Cable, User Manual

Product Description

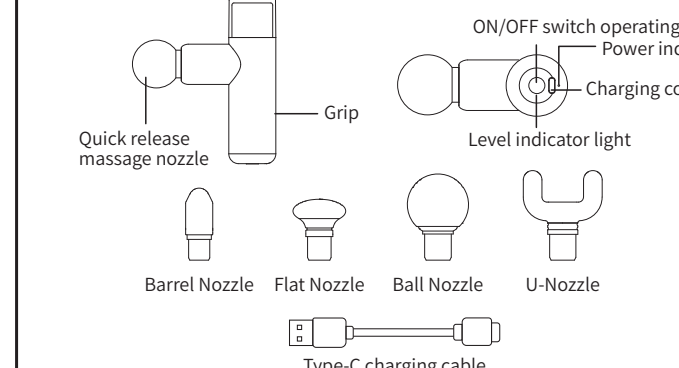
Product Operating Instructions

To charge this product, connect the USB-C charging cable to it, and then plug the USB connector of the charging cable into a charger (please use a charger that supports the 5V=2A specification).
Reminder : Do not leave the site during charging.

1. The red battery indicator light stays on during charging.
2. The green light stays on when fully charged.
3. This product does not support usage while charging.
4. During product usage, the green battery indicator light remains on.

How To Use

1. Before installing the massage head, please make sure to turn off the machine. Choose the appropriate massage head according to the different parts of the body, and plug in the product.

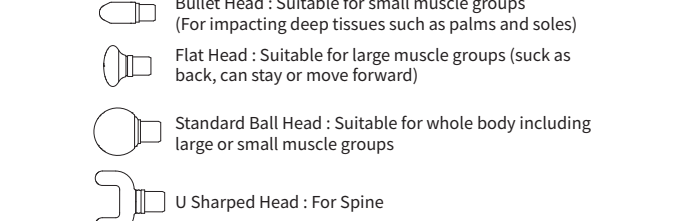2. Long press ON/OFF button for 2 sec. to power on. (Enter level 1 after power on)
3. Short press ON/OFF button to adjust level.

Level	Level indicator light	RPM (Revolutions per minute)
Level 1	The green light around the button is on	1600rpm
Level 2	The blue light around the button is on	2200rpm
Level 3	Blue and green mixed color lights are on	3000rpm
Level 4	Blue and green lights flash alternately	700rpm

Reminder : It is advised to use level 1 for beginner
4. When excessive pressure is applied to the massage gun, it will automatically activate a protection mechanism and pause operation. To resume use, simply briefly press the power button within 5 sec. If the pause exceeds 5 sec., the device will shut down automatically and require a restart to activate.
5. The product has smart timer function. It will stop working after 10 min. (Press ON/OFF button to reuse the product)
6. Long press ON/OFF button for 2 sec. to power off.

CAUTION

1. Do not use for prolonged periods. (Avoid continuous use for more than 15 sec at the same point; do not use on the same area for more than 10 min.)
2. Consult a doctor before using this device if you :
 - Need rest or have obvious physical discomfort or abnormal sensations.
 - Have peripheral circulation disorders such as diabetes-related sensory impairment.
 - Have a body temperature of 38°C (100.4°F) or higher.
 - Have wounds in the massage area.
 - Are diagnosed with malignant tumors.
 - Have heart conditions.
 - Have lost the ability to sense heat.
3. Avoid vibrating over bones and joints.

WARNING

1. This appliance is not suitable for unsupervised children or individuals with physical/mental disabilities.
2. Supervise children to prevent them from playing with the appliance.
3. Only use the appliance with the indicated extra-low safety voltage.
4. Do not use this product on areas of high vascular concentration on the body, including the head, face, and neck vertebrae.
5. Do Not Use if:
 - Advised against massage by a physician (e.g., thrombosis, severe aneurysm, acute varicose veins, skin disorders/infections).
 - Acute lower back pain, abnormal spine conditions, herniated discs, spinal misalignment, stenosis, or degenerative lumbar facet syndrome.
 - Osteoporosis, fractures, sprains, or acute pain from muscle strains/tears.
 - Implanted medical devices sensitive to electromagnetic interference (e.g., pacemakers).
 - Pregnant or recently given birth.

Safety Instruction For Electronic Product

1. During use, do not insert fingers, any part of the body, or other objects into the gap between the massage head and the product body to avoid pinching.
2. This product is not a substitute for medical treatment. If you have pain or are undergoing treatment, consult a physician before use.
3. This product is not waterproof. Do not use it in damp environments or allow the product to come into contact with or be submerged in water.
4. Clean with a soft cloth. Do not use strong chemical products or detergents to clean the device.
5. Avoid dropping, strong vibrations, or impacting the product with hard objects to prevent damage to internal circuits.
6. Do not place the product in direct sunlight or in high-temperature environments.
7. For safety, do not leave the product unattended while charging.
8. Do not use this product near any heat source.
9. Do not disassemble screws or dismantle the product by yourself.
10. Do not charge and use the product at the same time.
11. Do not replace the massage head while the product is in use.

Environmental Friendliness

Cleaning

1. Turn off and unplug.
2. Wipe exterior with a damp cloth.
3. Clean with soapy water or alcohol, then dry.
4. Air dry in a ventilated area.

Maintenance


1. Use for no more than 15 minutes at a time.
2. Store in a cool, dry place.
3. Follow charging recommendations.

Environmental Friendliness

This product has built-in rechargeable lithium battery. When discarding this product, do not dispose it together with household garbage. Be sure to safely handle this product at hazardous waste collection sites, recovery facilities or recycling electronics retailers.

Battery Disposal Instructions

Before discarding this product, the battery must be removed and disposed of safely. Disconnect the product before removing the battery. Remove the battery from the bottom, as shown in the diagram :

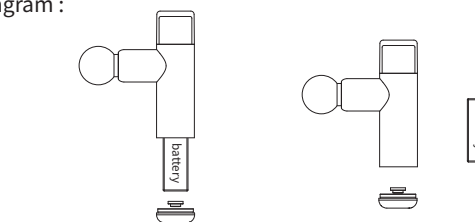Warranty Instructions

Warranty Period Description :
For VIP gifts such as mobile phones/tablets sold by Samsung Electronics Taiwan, you can enjoy a 1-year warranty from the date of purchase with this warranty and the purchase invoice of Samsung products.
The free service will not be provided in the following cases:

1. The product has exceeded the warranty period provided by the company.
2. Damage or failure caused by natural disasters, insect (rat) damage, lightning strikes and other factors.
3. Failures caused by other products (equipment) other than the machine itself.
4. Failure to operate in accordance with product specifications and user manuals, resulting in product abnormalities or operating system changes.

5. The cause of the failure is beyond the scope of the original warranty, such as : damp, damaged, deformed, using a non-original authorized repair center for the product or using non-original spare parts.
6. Human failure during movement, transportation or improper placement.
7. Failures caused by using other equipment and accessories that do not meet the applicable specifications of the product.
8. Unsuitable power circuit and voltage and current.
9. Appearance and consumable materials.

SAMSUNG C&T TAIWAN CORPORATION
Address: 110 Rim. 2203 22F, No.333, Keelung Rd., Sec.1 Taipei, Taiwan
TEL: +886 2 2728 8500
E-mail: itfittw@samsung.com